



Precision in Healthcare

HOLLOW FIBER HEMODIALYZER

透析器

BASIC PARAMETERS OF THE FILTER

Model	QB (ml/min)	The Clearance Rate (ml/min)				Ultrafiltration Rate (ml/h)	Blood volume (ml)	Effective membrane area (m ²)
		Urea	Creatinine	Phosphate	Vitamin B12			
HD-10H	200/300/400	175/212/236	135/188/202	130/152/162	78/83/99	400	70	1.0
HD-12H	200/300/400	180/227/255	166/206/228	141/165/188	88/92/112	450	78	1.2
HD-13H	200/230/400	183/234/266	171/212/236	145/170/200	95/108/120	500	79	1.3
HD-14H	200/230/400	186/238/270	175/215/241	148/175/205	97/110/121	550	85	1.4
HD-15H	200/250/400	186/241/278	176/221/248	150/184/220	99/112/125	600	87	1.5
HD-16H	200/300/400	188/245/278	178/224/251	155/190/222	100/115/130	600	91	1.6
HD-17H	200/390/400	189/248/284	180/229/258	162/196/230	100/116/130	800	100	1.7
HD-18H	190/250/280	180/250/280	181/230/260	164/195/232	110/130/140	800	103	1.8
HD-19H	191/253/290	191/253/290	184/236/268	164/200/246	115/132/145	850	109	1.9
HD-20H	192/255/290	192/255/290	187/236/268	175/202/247	118/134/145	850	113	2.0
HD-21H	193/256/290	193/256/290	188/238/270	176/204/248	120/136/146	900	117	2.1
HD-22H	194/257/290	194/257/290	188/240/272	176/206/250	121/137/147	950	119	2.2

Test condition: QB-500mL/min, QF-10mL/min, QB-300mL/min

SCREENING COEFFICIENT (S.C.)

Test conditions:

QB = 200mL/min

QF = 30mL/min

Inulin: 0.9

β 2-microglobulin: 0.8

Myoglobin: 0.223

Albumin: < 0.002

In clinical application, the results
maybe slightly different due to the
use of different ultrafiltration methods.


STABLE PHYSICAL AND CHEMICAL PROPERTIES

We use polyether sulfone membrane imported from Germany as dialyzer membrane, which has excellent biocompatibility, reduces the use of anticoagulant, and has a three-layer asymmetric structure, which effectively intercepts endotoxins and makes dialysis treatment safer and more comfortable.



MORE STABLE PHYSICAL PROPERTIES

Better molecular structure
with higher stability.

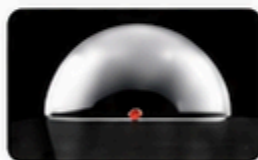


BETTER HYDROPHILICITY

Contact angle with water
is 67° (PES), 81° (PVP).



HIGHER THERMAL STABILITY AND CORROSION RESISTANCE

Glass transition temperature (T_g):
PES 230°C, PVP 190°C.

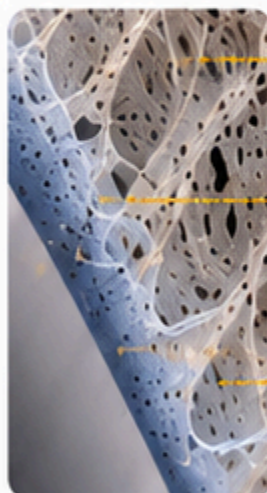

Extremely thin
endothelial layer
Ensures high dispersion
efficiency

Middle layer supports
the endothelium and
resists the entry of
endotoxins

Finger-shaped layer
ensures the strength
of the endothelium and
increases the efficiency of
diffusion and convective
transport



PES
MEMBRANE



S.E.T.
TECHNOLOGY



A.S.M.
TECHNOLOGY



P.E.T.
TECHNOLOGY



PYROGEN
FREE



STERILE R
EO / IR
STERILE



MANUFACTURER BY

JIANGXI HONGDA MEDICAL EQUIPMENT GROUP LTD

Shanghai International Holding Corp. GmbH(Europe)Hamburg Germany



MIPRO HEALTHCARE CORPORATION

Marketing by MIPRO HEALTHCARE CORPORATION



ISO 13485
CERTIFIED

CE
0123

STERILE R
EO / IR



PYROGEN
FREE

Trusted Technology
Better Life